



Package and Mail Inspection Workstation

- Provides Added Protection from Vapors and Particulates During Package Inspection




54 watt¹ Purair P5-36-XT (SAFESEARCH) shown with optional air velometer.



30–71 watt¹

The single EC blower motor assures lower cost of ownership in one of the world's most energy efficient ductless fume hoods.

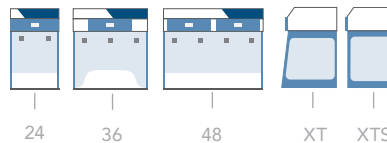

CONTENTS:

- Product Overview (p.2)
- Design Features (p.3)
- Performance & Selection (p.4)
- Filtration Technology (p.5)
- Specifications (p.6)
- Options & Accessories (p.8)

INTRODUCTION

Purair® SAFESEARCH ductless fume hoods are a series of high efficiency workstations designed to increase protection from drugs, chemical vapors and particulates during the packaging examination process. The innovative [Multiplex™ Filtration Technology](#) combines a broad spectrum, general purpose carbon filter with a HEPA filter to protect the investigator and the environment.

APPLICATIONS

Using innovative filtration technology, the Purair SAFESEARCH creates a safe work environment over the widest range of applications in the industry.

Routine or Suspicious Envelope and Package Inspection \ Forensics \ Evidence Examination \ Mailrooms \ Airports \ Bus Terminals \ Vendor Receiving Bays \ Government Offices \ Law Enforcement \ Correctional Facilities



Deep into its second generation, Air Science embraces the diversity and cultural heritage of the founders and co-workers who are continuing a tradition of excellence. Demonstrating a commitment to adaptation, inclusion and quality output from a United States-based company with a domestic and global reach.

Air Science® USA LLC \ 120 6th Street, Fort Myers, FL 33907
T. 239-489-0024 \ airscience.com

Air Science® Technologies Ltd. \ United Kingdom
T. 0151 526 2457 \ airscience.com/UK



KEY FEATURES

- Protects the investigator from hazardous drugs, chemicals and particulates.
- Easy to change filters.
- High efficiency EC blower.
- Improved filter clamping eliminates bypass leakage.
- Filter blockage alarm.

DUCTLESS TECHNOLOGY

The Eco-Friendly Choice

Advanced carbon filtration technology offers a safe, high performance alternative to conventional ducted fume hoods for a broad range of applications.

Environmental Benefits. Air Science® ductless fume hoods isolate and trap chemical vapors to prevent ecological impact through release into the environment.

Versatile. Each workstation includes a broad spectrum multi-chemical carbon filter combined with a HEPA filter. A range of carbon filter chemical families can be specified depending on the type of contamination that may be expected in a particular facility.

Easy to Install. The ductless fume hood is self-contained and does not require venting to the outside. Many units are portable and may be moved from one location to the next with minimal downtime and without filter changes. Set-up, operation and filter maintenance are straightforward.

Energy Efficient. Because filtered air is returned to the room, no demands are required of the facility HVAC capacity for make-up air.

Cost Effective. Facility ductwork, HVAC and construction costs are eliminated.

Safe to Use. Cabinet airflow and face velocity protect users from incidental exposures to fumes.

Self-Testing. (selected models) Electronic airflow monitoring ensures continuous safety. An electronic gas sensor monitors carbon filter performance.

Specifications are subject to change without notice or obligation on the part of Air Science. For questions contact Air Science.

¹⁾ Energy consumption disclosure is based on internal testing with primary filters during normal operation. Power consumption published is nominal and dependent on cabinet size.



71 watt¹ Purair P5-48-XT (SAFESEARCH) shown with optional velometer and mobile base stand with casters.



30–71 watt¹

The single EC blower motor assures lower cost of ownership in one of the world's most energy efficient ductless fume hoods.

CONTENTS:

- Product Overview (p.2)
- Design Features (p.3)
- Performance & Selection (p.4)
- Filtration Technology (p.5)
- Specifications (p.6)
- Options & Accessories (p.8)

PurairSAFESEARCH

Ductless Fume Hoods

24 • 36 • 48

DESIGN FEATURES

3



DESIGN FEATURES

- A. Filter I.D. Window:** A convenient, strategically placed front cover window shows the installed filter part number and installation date to encourage timely filter replacement.
- B. Control Panel:** Electronic controls and displays include switches for the blower and filter blockage alarm.
- C. Filter Blockage Alarm:** Continuously monitors filter loading and alerts the user when service is needed.
- D. Air Velometer:** An optional analog air velocity meter is positioned in the user's field of vision.
- E. Steel Support Frame:** The chemical resistant epoxy coated steel frame adds mechanical strength.
- F. Hinged Front Sash:** The sash is easy to open and close. When closed, the cabinet sash protects the contents from inadvertent external contact and better isolates the air within.
- G. Work Surface:** The internal work surface includes a black, polypropylene work tray to improve visibility of white powders such as illegal drugs.
- H. Pass Through Ports:** Electrical cords and cables are safely routed into the cabinet through ports on the back.
- I. Electrostatic Pre-Filter:** The 99.5% effective electrostatic pre-filter is accessible from inside the chamber to contain the release of any particulates that it traps. The pre-filter can be changed while the unit is operating to prevent operator exposure to chemical vapors.

- J. Filter Door Key:** Filter access keys prevent unauthorized removal or accidental exposure to dirty filters.
- K. Internal Manual Speed Controller:** Authorized personnel may set the EC blower speed as desired.
- L. Stand:** Optional mobile base stand with locking casters.

ADDITIONAL FEATURES

360 Degree Visibility: Clear back and side panels allow ambient light into the chamber and provide users with an unobstructed view of contents.

Standards Compliant: Performance specifications and construction meet or exceed relevant standards to ensure operator safety.



30 watt¹ Purair P5-24-XT (SAFESEARCH), shown with optional velometer and mobile base stand with casters. Ductless filtration permits easy relocation from one mailroom or processing facility to another.

Air Science® USA LLC \ 120 6th Street, Fort Myers, FL 33907
T. 239-489-0024 \ airscience.com

Air Science® Technologies Ltd. \ United Kingdom
T. 0151 526 2457 \ airscience.com/UK

Specifications are subject to change without notice or obligation on the part of Air Science. For questions contact Air Science.

¹) Energy consumption disclosure is based on internal testing with primary filters during normal operation. Power consumption published is nominal and dependent on cabinet size.

CONTENTS:

- Product Overview (p.2)
- Design Features (p.3)
- Performance & Selection (p.4)
- Filtration Technology (p.5)
- Specifications (p.6)
- Options & Accessories (p.8)

PurairSAFESEARCH

Ductless Fume Hoods

24 • 36 • 48

PERFORMANCE & SELECTION

4

Air Science high-efficiency fume hoods are expertly designed to meet specific applications and certified for quality construction. Standard features, options and accessories are developed purposefully to enhance user-friendliness.

PERFORMANCE

The Purair SAFESEARCH accommodates the full range of [Multiplex Filtration System](#) options.

The high capacity air handling system delivers face velocity of 100 fpm in compliance with US and international safety and performance standards.

DESIGN

Professional quality Air Science fume hoods comply with current technical and safety regulations.

The cabinet frame and work surfaces, comprised of industrial components, are durable and chemically resistant.


The Air Science filter assembly is easy to access and change.

RELIABILITY

Internal systems are isolated from fumes, extending product life.



Energy-efficient blowers promote long life and dependable performance of Purair SAFESEARCH fume hoods.

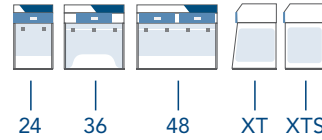

SELECTION

Purair SAFESEARCH products are available in 3 standard widths and 2 depth options, in metal construction, totaling 6 standard models.


CONTROL

The **Basic control panel** is standard on Purair SAFESEARCH models and includes an On/Off switch and filter blockage alarm.

The **optional Advanced controller** displays the airflow and offers limited detection of low concentrations of hydrocarbon, some gases and organic acids. Audio and visual alarms alert users to filter saturation and if the airflow reaches preset thresholds. An Hour Counter is also included.



Basic Control Panel

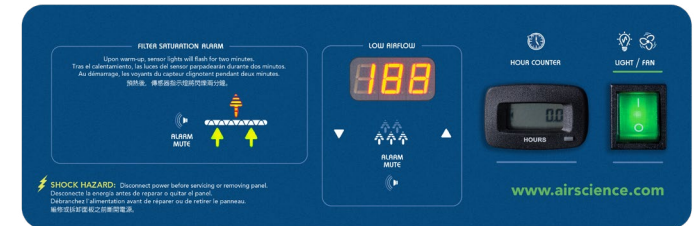

Advanced Control Panel (ADVP)

Air Science® USA LLC \ 120 6th Street, Fort Myers, FL 33907
T. 239-489-0024 \ airscience.com

Air Science® Technologies Ltd. \ United Kingdom
T. 0151 526 2457 \ airscience.com/UK

Specifications are subject to change without notice or obligation on the part of Air Science. For questions contact Air Science.

¹⁾ Energy consumption disclosure is based on internal testing with primary filters during normal operation. Power consumption published is nominal and dependent on cabinet size.

CONTENTS:

- Product Overview (p.2)
- Design Features (p.3)
- Performance & Selection (p.4)
- Filtration Technology (p.5)**
- Specifications (p.6)
- Options & Accessories (p.8)



FILTRATION

At the heart of the Purair product line is innovative filtration technology. **The Multiplex Filtration System** consists of a pre-filter and main carbon and HEPA filters. The system permits a customized combination of filter media and configuration for chemical and physical adsorption specific to each application need.

The carbon filter media is treated to attain the proper porosity and aggregate surface area and to react with several ranges of aerosolized chemicals while the integrated HEPA filter captures particulates, powders, or biological contaminants. The mechanical design enhances safety, convenience and overall value.

View available filters and descriptions on [page 7](#).



Filter disposal services are available in selected markets providing responsible destruction or recycling of used saturated filters in authorized facilities.

FILTER CONFIGURATION

The Multiplex feature permits one or more filtration options to be combined to meet a wider range of multiple-use applications. The Purair SAFESEARCH can be equipped with a single activated carbon main filter or with a stacked configuration which combines two main filters.

The carbon filter is sized to fit the specified product model number and configured to optimize airflow across 100% of the filter surface area. The self-contained assembly maximizes filter efficiency, prolongs filter life, optimizes diffusion and saturation and improves user safety.

- P. Electrostatic Pre-Filter:** Protects the main filters from aerosols, mists, dust and particulates.
- C. Activated Carbon Main Filter:** A single or stacked filter configuration.
- H. HEPA/ULPA Filter:** Both HEPA and optional ULPA filters use micro-glass fiber media designed to capture fine particles and biologicals. Both filters can capture particles smaller than the micron size for which they are tested. HEPA and ULPA filter efficiencies are 99.97% at 0.3 microns and 99.999% at 0.12 microns respectively.

MULTIPLY FILTRATION SYSTEM, SUMMARY				
Application	Chemical	Powder/ Biological	Chemical & Powder	Chemical within Cleanroom
Secondary/ Stacked Filter, Optional	C	H	C	H
Primary Filter	C	H	H	C
Pre-Filter	P	P	P	P


The system can be configured for the capture of acids, bases and particulates, such as biological aerosols, when paired with HEPA or ULPA filters.

AIRFLOW

Contaminated air is pulled through the Multiplex Filtration System where activated carbon adsorbs chemical vapors and the HEPA filter captures particulates. Clean air is returned to the room.

The main filters are easy to replace, no tools required. The filter clamps tightly against the filter gasket to prevent filter bypass and maintain filter integrity.

⚠ The pre-filter may be replaced while unit is in operation.

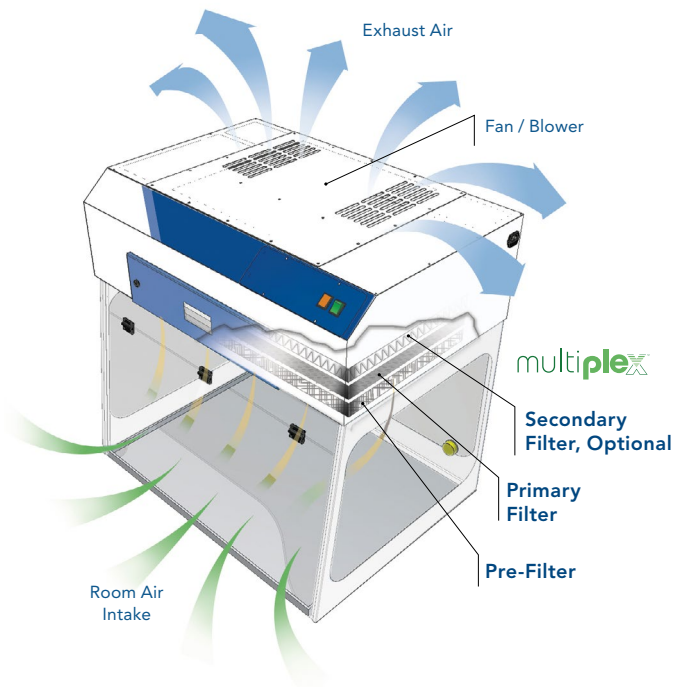


Specifications are subject to change without notice or obligation on the part of Air Science. For questions contact Air Science.

¹⁾ Energy consumption disclosure is based on internal testing with primary filters during normal operation. Power consumption published is nominal and dependent on cabinet size.


CONTENTS:

- Product Overview (p.2)
- Design Features (p.3)
- Performance & Selection (p.4)
- Filtration Technology (p.5)
- Specifications (p.6)
- Options & Accessories (p.8)


**P5-24-XT
(SAFESEARCH)**

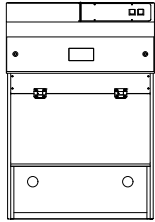


**P5-36-XT
(SAFESEARCH)**

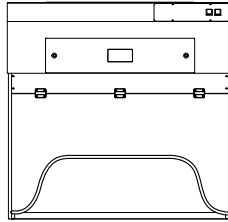


**P5-48-XT
(SAFESEARCH)**

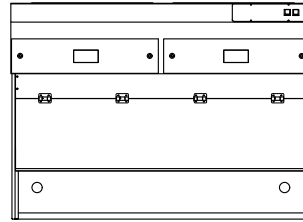

Standard Depth (XT)

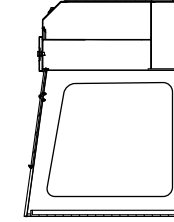

Reduced Depth (XTS)

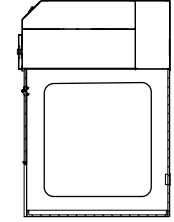

MODEL	VOLTAGE	DIMENSIONS			WEIGHT (LBS/KG)	
		Internal Height	External (W × D × H)	Shipping (W × D × H)	Net	Ship
Standard Depth Models (XT)						
P5-24-XT-A (SAFESEARCH)	120V, 60Hz	23.6" / 600 mm	24" × 27" × 35" / 610 × 676 × 889 mm	40" × 40" × 40" / 1016 × 1016 × 1016 mm	72 / 33	150 / 68
P5-24-XT-G (SAFESEARCH)	230V, 50Hz	23.6" / 600 mm	24" × 27" × 35" / 610 × 676 × 889 mm	40" × 40" × 40" / 1016 × 1016 × 1016 mm	72 / 33	150 / 68
P5-36-XT-A (SAFESEARCH)	120V, 60Hz	23.6" / 600 mm	36" × 27" × 35" / 914 × 676 × 889 mm	40" × 40" × 40" / 1016 × 1016 × 1016 mm	99 / 45	170 / 77
P5-36-XT-G (SAFESEARCH)	230V, 50Hz	23.6" / 600 mm	36" × 27" × 35" / 914 × 676 × 889 mm	40" × 40" × 40" / 1016 × 1016 × 1016 mm	99 / 45	170 / 77
P5-48-XT-A (SAFESEARCH)	120V, 60Hz	23.6" / 600 mm	48" × 27" × 35" / 1219 × 676 × 889 mm	55" × 40" × 45" / 1397 × 1016 × 1143 mm	138 / 63	230 / 104
P5-48-XT-G (SAFESEARCH)	230V, 50Hz	23.6" / 600 mm	48" × 27" × 35" / 1219 × 676 × 889 mm	55" × 40" × 45" / 1397 × 1016 × 1143 mm	138 / 63	230 / 104
Reduced Depth Models (XTS) for countertops 24" or less						
P5-24-XTS-A (SAFESEARCH)	120V, 60Hz	23.6" / 600 mm	24" × 24" × 35" / 610 × 610 × 889 mm	40" × 40" × 40" / 1016 × 1016 × 1016 mm	72 / 33	150 / 68
P5-24-XTS-G (SAFESEARCH)	230V, 50Hz	23.6" / 600 mm	24" × 24" × 35" / 610 × 610 × 889 mm	40" × 40" × 40" / 1016 × 1016 × 1016 mm	72 / 33	150 / 68
P5-36-XTS-A (SAFESEARCH)	120V, 60Hz	23.6" / 600 mm	36" × 24" × 35" / 914 × 610 × 889 mm	40" × 40" × 40" / 1016 × 1016 × 1016 mm	99 / 45	170 / 77
P5-36-XTS-G (SAFESEARCH)	230V, 50Hz	23.6" / 600 mm	36" × 24" × 35" / 914 × 610 × 889 mm	40" × 40" × 40" / 1016 × 1016 × 1016 mm	99 / 45	170 / 77
P5-48-XTS-A (SAFESEARCH)	120V, 60Hz	23.6" / 600 mm	48" × 24" × 35" / 1219 × 610 × 889 mm	55" × 40" × 45" / 1397 × 1016 × 1143 mm	138 / 63	230 / 104
P5-48-XTS-G (SAFESEARCH)	230V, 50Hz	23.6" / 600 mm	48" × 24" × 35" / 1219 × 610 × 889 mm	55" × 40" × 45" / 1397 × 1016 × 1143 mm	138 / 63	230 / 104

CONTENTS:

- Product Overview (p.2)
- Design Features (p.3)
- Performance & Selection (p.4)
- Filtration Technology (p.5)
- Specifications (p.6)
- Options & Accessories (p.8)

PRODUCT SPECIFICATIONS

Filtration	P5-24-XT (SAFESEARCH) P5-24-XTS (SAFESEARCH)	P5-36-XT (SAFESEARCH) P5-36-XTS (SAFESEARCH)	P5-48-XT (SAFESEARCH) P5-48-XTS (SAFESEARCH)
Face Velocity	100 fpm	100 fpm	100 fpm
Construction	P5-24-XT (SAFESEARCH) P5-24-XTS (SAFESEARCH)	P5-36-XT (SAFESEARCH) P5-36-XTS (SAFESEARCH)	P5-48-XT (SAFESEARCH) P5-48-XTS (SAFESEARCH)
Finish	<... White epoxy coated steel frame and head unit. Clear sides and back panel. ...>		
Blower	<... EC blower. ...>		
Controls	<... Main On/Off. ...>		
Monitoring	<... Filter blockage alarm, standard. ...>		
Efficiency	P5-24-XT (SAFESEARCH) P5-24-XTS (SAFESEARCH)	P5-36-XT (SAFESEARCH) P5-36-XTS (SAFESEARCH)	P5-48-XT (SAFESEARCH) P5-48-XTS (SAFESEARCH)
Power Consumption* 120V, 60Hz	30 watt	54 watt	71 watt
Power Consumption* 230V, 50Hz	30 watt	54 watt	71 watt
Lighting	<... LED. ...>		

* Watts at calibrated airflow setpoint with GP Filter(s) and prefilter installed.

FILTER SPECIFICATIONS

Purair Model	P5-24-XT (SAFESEARCH) P5-24-XTS (SAFESEARCH)	P5-36-XT (SAFESEARCH) P5-36-XTS (SAFESEARCH)	P5-48-XT (SAFESEARCH) P5-48-XTS (SAFESEARCH)
Secondary/Stacked Filter, Optional*	(1)	(1)	(2)
Primary Filter*	(1)	(1)	(2)
Pre-Filter*	(1)	(1)	(2)

* For specific examples refer to Multiplex filtration system summary on [page 5](#).

FILTER SUMMARY

Formula	Description
GP Plus!	The most widely used filter in the range, primarily for solvent, organic and alcohol removal.
HEPA	Powders and particulates.

View additional information on the Multiplex Filtration System on [page 5](#).



Through our partner company [Filtco Filters](#), Air Science is a single source supplier of all pre-filters, carbon filters and HEPA/ULPA filters used in our products.

CONTENTS:

- Product Overview (p.2)
- Design Features (p.3)
- Performance & Selection (p.4)
- Filtration Technology (p.5)
- Specifications (p.6)
- Options & Accessories (p.8)

OPTIONS & ACCESSORIES

Purair Model		P5-24-XT (SAFESEARCH) P5-24-XTS (SAFESEARCH)	P5-36-XT (SAFESEARCH) P5-36-XTS (SAFESEARCH)	P5-48-XT (SAFESEARCH) P5-48-XTS (SAFESEARCH)
Advanced Controller*	The optional Advanced controller displays the airflow and offers limited detection of low concentrations of hydrocarbon, some gases and organic acids. Audio and visual alarms alert users to filter saturation and if the airflow reaches preset thresholds. An Hour Counter is also included.	P5-24-XT-A-ADVP (SAFESEARCH) P5-24-XT-G-ADVP (SAFESEARCH) P5-24-XTS-A-ADVP (SAFESEARCH) P5-24-XTS-G-ADVP (SAFESEARCH)	P5-36-XT-A-ADVP (SAFESEARCH) P5-36-XT-G-ADVP (SAFESEARCH) P5-36-XTS-A-ADVP (SAFESEARCH) P5-36-XTS-G-ADVP (SAFESEARCH)	P5-48-XT-A-ADVP (SAFESEARCH) P5-48-XT-G-ADVP (SAFESEARCH) P5-48-XTS-A-ADVP (SAFESEARCH) P5-48-XTS-G-ADVP (SAFESEARCH)
Dwyer Airflow Meter	Continuous display of face velocity.	DWYER	DWYER	DWYER
Base Stand, Mobile, with Casters	Mobile base stand, fixed height, with locking casters.	CART-25	CART-36	CART-50
Bottom Shelf	Provides a lower storage shelf for mobile base stand.	CART-25-SHELF	CART-36-SHELF	CART-50-SHELF
Base Cabinet, Fixed (Metal)	Provides storage space below.	CART-MCC-25	CART-MCC-36	CART-MCC-50
Base Cabinet, Fixed (Polypropylene)	Provides storage space below.	CART-SSC-25	CART-SSC-36	CART-SSC-50
Fire Safety Cabinet Base	Flame resistant safe storage for combustible and flammable liquids.	CART-FSC-25	CART-FSC-36	CART-FSC-50

*Specify when ordering.

"A" — 120V, 60Hz

"G" — 230V, 50Hz

Certain options, customizations or configurations may not be included in UL-C-61010-1 listings. Contact Air Science for details. Specifications are subject to change without notice or obligation on the part of Air Science. For questions contact Air Science.

CONTENTS:

- Product Overview (p.2)
- Design Features (p.3)
- Performance & Selection (p.4)
- Filtration Technology (p.5)
- Specifications (p.6)
- Options & Accessories (p.8)

WARRANTY

This product is protected by the Air Science Legacy Limited Lifetime Warranty™.



For details visit the [Warranty section](#) of our website.

STANDARDS & COMPLIANCE

Quality Management Systems	ISO 9001: 2015
Electrical Safety	UL-C-61010-1 CAN/CSA C22.2 61010-1-12 EN 61010-1:2010 CE Mark
OSHA, Occupational Safety and Health Administration	OSHA Standard -29 CRF, Safety and Health Regulations for General Industry, 1910.1450: Occupational exposure to hazardous chemicals in laboratories. Part B, definition, laboratory type hood. This product may assist you with compliance or as part of your chemical hygiene plan. Please consult your Safety Officer and/or Industrial Hygienist.
Environment	ISO 14001: 2015 ENERGY STAR® Partner



120 6th Street \ Fort Myers, FL 33907
T. 239-489-0024 \ Toll Free. 800-306-0656 \ F. 800-306-0677
airscience.com

Air Science® Technologies Ltd. \ United Kingdom
T. 0151 526 2457 \ airscience.com/UK

The information contained in this manual and the accompanying product are copyrighted and all rights are reserved by Air Science. Air Science reserves the right to make periodic minor design changes without obligation to notify any person or entity of such change.

